


XA 4Y Analytical Balances

Analytical weighing of the highest quality and accuracy using professional balances





XA 4Y, d = 0.1 mg

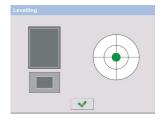

Automatic control of the level



Intuitive operation and touch

Functions



counting

Dosing

Checkweighing

Formulations

Percent weighing

Statistics

weighing

Differential weighing

Pipettes calibration

Statistical quality control

Autotest



Density determination

Under hook weighing

Cooperation with titrators

Ambient conditions measurement



GLP procedures



Proximity sensors



Moveable range



Replaceable unit



Multilingual

Features

Reliable Results and Excellent Weighing Parameters

The best weighing parameters enable applying XA 4Y analytical balances in majority of demanding laboratory processes.

Innovative 2-Point Adjustment

The unique 2-point adjustment system guarantees the best accuracy of the weighing results and minimizes linearity errors thus ensuring reliable measurements within the whole measuring range.

Significantly Fast Measurements

Powerful processor offers new possibilities of operation assuring short indication stabilization time at a respective repeatability.

Second to None Repeatability and Compliance with USP

XA 4Y analytical balances feature the highest measurements accuracy, excellent repeatability and are compliant with USP requirements (Chapter 41 and 1251).

Intuitive Operation and Touch Screen

5.7" colour touch screen enables intuitive operation and easy access to numerous applications and functions of the weighing instrument.

Automatic Control of the Level

Levelling system facilitates adjustment of device level, it also uninterruptedly controls the level state, and informs about potential level deviations.

Spacious Weighing Chamber

Large weighing chamber enables convenient operation using laboratory vessels of different dimensions.

Numerous Options of Data Management

Extensive storage capacity allows record of all measurement data in a form of complex reports and statistical graphs.

Page 1 of 5 | Date: 21.08.2018 www.radwag.com

Technical Specifications

	XA 52.4Y	XA 82/220.4Y	XA 110.4Y	XA 120/250.4Y
Maximum capacity [Max]	52 g	82/220 g	110 g	120 g /250 g
Minimum load	1 mg	1 mg	1 mg	1 mg
Readability [d]	0.01 mg	0.01 mg / 0.1 mg	0.01 mg	0.01g /0.1 mg
Verification scale interval [e]	1 mg	1 mg	1 mg	1 mg
Tare range	-52 g	–220 g	–110 g	-250 g
Repeatability(5% Max)*	0.006 mg (Rt ≤ 3g)	0.007 mg (Rt ≤ 10 g)	0.007 mg (Rt ≤ 5 g)	0.007 mg (Rt ≤ 10 g)
Repeatability(Max)*	0.01 mg	0.06 mg	0.02 mg	0.06 mg
Linearity	±0.03 mg	±0.06 / ±0.2 mg	±0.06 mg	±0.06 / ±0.2 mg
Eccentric load deviation	0.03 mg	0.2 mg	0.06 mg	0.2 mg
Sensitivity temperature drift**	1×10^{-6} / °C × Rt	1×10^{-6} /°C × Rt	1 × 10 ⁻⁶ / °C × Rt	1×10^{-6} /°C×Rt
Sensitivity time drift	1×10^{-6} / Year \times Rt			
Minimum weight (U=1%, k=2)	1.2 mg	1.4 mg	1.4 mg	1.4 mg
Minimum weight (USP)	12 mg	14 mg	14 mg	14 mg
Stabilization time	~ 4 s	~ 4 s	~ 4 s	~ 4 s
Adjustment	internal	internal	internal	internal
Moveable range	_	Yes	_	Yes
Verification	Yes	Yes	Yes	Yes
OIML Class	1	1	1	1
Indicator fastening	35 cm cable, wireless connection (option)***			
Display	5.7" colour, resistive touch screen			
Keypad	8 keys	8 keys	8 keys	8 keys
Protection class	IP 43	IP 43	IP 43	IP 43
Databases	19	19	19	19
Touch-free operation	2 programmable proximity sensors			
USB-A	2	2	2	2
Ethernet	10 / 100 Mbit			
RS 232	2	2	2	2
Wireless connection	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n
IN/OUT	$4 \times IN, 4 \times OUT$			
Power supply	13.5 ÷ 16 V DC			
Power consumption	10 W	10 W	10 W	10 W
Operating temperature	+10 ÷ +40 °C			
Atmospheric humidity****	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%
Transport and storage temperature	-20 ÷ +50 °C	-20 ÷ +50 °C	−20 ÷ +50 °C	-20 ÷ +50 °C
Weighing pan dimensions	ø 90 mm open-work ø 85 mm standard (option)*****	ø 90 mm open-work ø 85 mm standard (option)*****	ø 90 mm open-work ø 85 mm standard (option)*****	ø 90 mm open-work ø 85 mm standard (option)*****
Weighing chamber dimensions	168 × 160 × 228 mm			
Weighing device dimensions	542 × 206 × 303 mm			
Net weight	9.8 kg	9.8 kg	9.8 kg	9.8 kg
Gross weight	14.3 kg	14.3 kg	14.3 kg	14.3 kg
Packaging dimensions	720 × 385 × 485 mm			

Rt net weight

Values of parameters provided in Technical Specifications table have been determined under stable laboratory conditions. Due to ambient conditions impact or/and balance setup, the above parameters may vary for environments other than laboratory.

www.radwag.com

^{*} repeatability is expressed as a standard deviation from 10 weighing cycles

 $^{^{**}}$ parameter determined in the following temperature range: $+15 \div +35$ °C

^{***} optional solution on purchase order

^{****} non-condensing conditions

^{*****} ø 85 mm standard weighing pan on purchase order

	XA 210.4Y	XA 220.4Y	XA 310.4Y
Maximum capacity [Max]	210 g	220 g	310 g
Minimum load	1 mg	10 mg	10 mg
Readability [d]	0.01 mg	0.1 mg	0.1 mg
Verification scale interval [e]	1 mg	1 mg	1 mg
Tare range	–210 g	–220 g	–310 g
Repeatability(5% Max)*	0.007 mg (Rt ≤ 10g)	0.06 mg (Rt ≤ 10g)	0.07 mg (Rt ≤ 15g)
Repeatability(Max)*	0.03 mg	0.1 mg	0.15 mg
Linearity	±0.1 mg	±0.2 mg	±0.3 mg
Eccentric load deviation	0.1 mg	0.2 mg	0.3 mg
Sensitivity temperature drift**	1 × 10 ⁻⁶ / °C × Rt	1 × 10 ⁻⁶ / °C × Rt	1 × 10 ⁻⁶ / °C × Rt
Sensitivity time drift	1×10^{-6} / Year \times Rt	1×10^{-6} / Year \times Rt	1×10^{-6} / Year \times Rt
Minimum weight (U=1%, k=2)	1,4 mg	12 mg	14 mg
Minimum weight (USP)	14 mg	120 mg	140 mg
Stabilization time	~ 4 s	~ 2.5 s	~ 2.5 s
Adjustment	internal	internal	internal
Moveable range	_	_	_
Verification	Yes	Yes	Yes
OIML Class	T	I	I
Indicator fastening	35 cm cable. wireless connection (option)***	35 cm cable. wireless connection (option)***	35 cm cable. wireless connection (option)***
Display	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen	5.7" colour, resistive touch screen
Keypad	8 keys	8 keys	8 keys
Protection class	IP 43	IP 43	IP 43
Databases	19	19	19
Touch-free operation	2 programmable proximity sensors	2 programmable proximity sensors	2 programmable proximity sensors
USB-A	2	2	2
Ethernet	10 / 100 Mbit	10 / 100 Mbit	10 / 100 Mbit
RS 232	2	2	2
Wireless connection	802.11 b/g/n	802.11 b/g/n	802.11 b/g/n
IN/OUT	$4 \times IN. 4 \times OUT$	$4 \times IN.4 \times OUT$	$4 \times IN.4 \times OUT$
Power supply	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC	13.5 ÷ 16 V DC
Power consumption	10 W	10 W	10 W
Operating temperature	+10 ÷ +40 ℃	+10 ÷ +40 °C	+10 ÷ +40 °C
Atmospheric humidity****	40 ÷ 80%	40 ÷ 80%	40 ÷ 80%
Transport and storage temperature	−20 ÷ +50 °C	−20 ÷ +50 °C	−20 ÷ +50 °C
Weighing pan dimensions	ø 90 mm open-work ø 85 mm standard (option)*****	ø 100 mm standard	ø 100 mm standard
Weighing chamber dimensions	168 × 160 × 228 mm	168 × 160 × 228 mm	168 × 160 × 228 mm
Weighing device dimensions	542 × 206 × 303 mm	542 × 206 × 303 mm	542 × 206 × 303 mm
Net weight	9.8 kg	9.8 kg	9.8 kg
Gross weight	14.3 kg	14.3 kg	14.3 kg
Packaging dimensions	720 × 385 × 485 mm	720 × 385 × 485 mm	720 × 385 × 485 mm

***** ø 85 mm standard weighing pan on purchase order


Values of parameters provided in Technical Specifications table have been determined under stable laboratory conditions. Due to ambient conditions impact or/and balance setup, the above parameters may vary for environments other than laboratory.


www.radwag.com Page 3 of 5 | Date: 21.08.2018


repeatability is expressed as a standard deviation from 10 weighing cycles parameter determined in the following temperature range: +15 \div +35 $^{\circ}$ C

optional solution on purchase order ***


non-condensing conditions

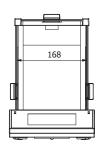


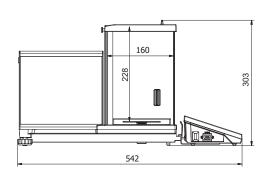


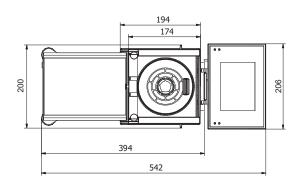



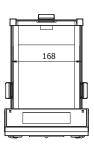

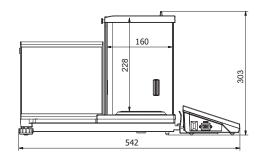

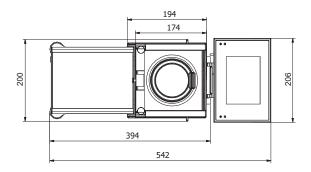

XA 4Y, d = 0.1 mg

Accessories

Weighing Tables

- granite antivibration table
- antivibration tables for laboratory balances
- professional weighing table

Professional Weighing

- Adapter for calibration of XA17 i XA100 series pipettes
- KIT 85 density determination kit
- · laboratory ware holders
- automatic feeder
- under-hook weighing rack

Ambient Conditions

- DJ-04 anti-static ioniser
- •THB-S or THB-P sensor

Peripheral Devices

- Epson dot matrix printer
- barcode scanners
- WD-5/3Y LCD display

Cables, Converters

- P0108: RS 232 cable (balance-computer)
- P0167: RS 232 cable (balance-computer)
- P0151: RS 232 cable (balance Epson printer)

Electrical Accessories

• ZR-02 power supply with battery

Dedicated Software

R-LAB

- collecting measurements
- carrying out statistical analysis of measurements
- customized graphs and reports

E2R Weighing Records

- complete, automated databases synchronization
- fully supported processes of labelling and parts counting
- record of weighings, weighings archiving
- · basic and advanced (with graphs) reports

Label Editor R02

- designing label templates
- sending graphics and fonts to label printers
- printing label templates using connected printers

Pipettes

- determining measurement errors of pipettes volume
- accordance with ISO 8655
- calibration of single-channel and multi-channel pipettes
- calibration of fixed-volume and variable-volume pipettes

Audit Trail Reader

- support of Audit Trail function available for 3Y, 4Y, HY10, WLY, WPY series weighing instruments
- record of operator's activity from the moment of logging in

Parameters Editor

- remote change of parameters
- remote on-line preview of the display
- displaying current mass indication
- software update
- file loading, editing and saving parameters to a file
- import and export of parameters
- interfaces: RS232, Ethernet and Wireless Connection.
- quick and easy edition of balance parameters using computer.

RAD KEY

• Establishing cooperation between a weighing instrument and a computer

R. Barcode

• The basic function software is presentation of the data sent by barcode scanners connected to PC via USB or RS232

Radwag Development Studio

- presentation of functions (and subfunctions) of communication protocol (Common Communication Protocol)
- possibility of connection with weighing equipment on which each function is carried out,
- library with mass control, contained within the development environment
- complete documentation of the communication protocol
- set of user manuals for different solutions addressed for programmers employed in companies using RADWAG-manufactured weighing equipment

LabView Driver

• operation of RADWAG balances in LabView environment

RADWAG Connect

- establishing communication with all balances, scales and weighing modules using Common Communication Protocol
- · communication via local network,
- support of basic functions
- · auto searching for devices
- connecting with few devices simultaneously, swapping between them
- clear list of connected platforms
- record of measurements in the program,
- export of carried out measurements to CSV file,
- work performed using freely selected device with Windows 10 operating system

RADWAG Remote Desktop

- remote operation via computer, mobile phone or tablet
- sending text messages
- version for Windows 10 and Android systems

Page 5 of 5 | Date: 21.08.2018 www.radwag.com